

WHEELCHAIRS

PURETILT

The patented PureTilt Tilt-in-Space Wheelchair is designed by an Occupational Therapist to provide safe and independent transitions to facilitate optimal comfort for a patient.

A manual tilt wheelchair that has a user-adjusting tilt control so that the wheelchair user can readjust their own posture throughout the day. This is vital for comfort and pressure relief.

Features:

- Height-adjustable, ergonomic, contoured back with adjustable contoured headrest
- User / caregiver-operated seat-tilt adjustment controlled by a built-in lever that uses a gas spring piston
- Padded, hinged, non-detachable arms swing up and back to facilitate transfer
- Back folds forward to facilitate storage and transport
- 22-inch rear wheels with aluminum hand rims for easy self- propulsion
- Two front casters swivel for easy maneuvering
- Rear anti-tippers standard
- Angle-adjustable leg rest assembly with:
  - User-operated mechanism with telescoping adjustable length
  - Footplate folds flat and locks in place
- HCPCS code E1161
- 300 lb (136.08 kg) Weight Capacity EVENLY DISTRIBUTED



Accessories: sold separately

- Item No. PT3000-SP — Right & Left Side Pads
- Item No. PT3000-FP — Larger Leg Pad & Foot Pad






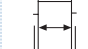



side pad set and foot pad accessories shown left  
- each sold separately



folds for storage and vehicle transport

Specifications:

Specifications:

Storage Dimensions with back folded										Seat Height front of seat to floor								Main, Seat, and Back Frame Wheels, Legrest, Footplate Other Durable Components <sup>4</sup>	
Overall Chair Width		Overall Chair Length				Length		Height with Seat		Seat Width		Seat Depth		bottom pin		top pin			
																			
Item	inch	cm	inch	cm	inch	cm	inch	cm	inch	cm	inch	cm	inch	cm	inch	cm	inch		cm
PT3000-18	28.50	72.39	39.80	101.09	50.00	127.00	39.50	100.33	36.00	91.44	19.80	50.29	21.00	53.34	19.75	50.17	21.75		55.25

Arm Height		Foot Support Length distance from floor				Front Caster Diameter		Rear Wheel Diameter		Seat Tilt adjustment <sup>2</sup>		Shipping <sup>3</sup> Weight		Product Weight with cushions		Weight Capacity EVENLY DISTRIBUTED		warranty (months)		
Item	inch	cm	retracted	cm	extended	cm	inch	cm	inch	cm	degree of movement	lb	kg	lb	kg	lb	kg	60	12	6
PT3000-18	5.80	14.73	13.00	33.02	17.50	44.45	6.00	15.20	22.00	55.90	5 to 40 -97 to 15	100	45.36	95	43.09	300	136.08	60	12	6

<sup>1</sup> ± 0.5 inches (± 1.27 cm) <sup>2</sup> ± 2° <sup>3</sup> estimated <sup>4</sup> not otherwise listed <sup>5</sup> from horizontal