

iM70

Veterinary Monitor



- 12.1-inch high-resolution color TFT-LCD screen
- Full touch screen
- Customizable shortcuts for simple and intuitive operation
- No-fan design reduces dust and creates a no noise environment
- Defibrillation and electrosurgical interference protection
- Basic parameters monitoring
- Optional capnography and anesthesia monitoring available
- Suitable for feline, canine and other animals

Standard Parameters: 3/5-lead ECG, HR, RESP, EDAN SpO₂, EDAN NIBP, PR, 2-TEMP

Optional Parameters: Nellcor OxiMaxTM SpO₂, SunTech NIBP, 2-IBP, EDAN G2 CO₂ (Sidestream), Respironics CO₂ (Mainstream and Sidestream), Masimo Anesthesia Gas/O₂ (iM70 VET)

1200

NIBP
Measurement

120 h

Trend
Review

60

Alarm
Review

120 s

Frozen
Waveform



Defib Proof & Sync



ESU Proof



VGA Output



iM70

12.1" Color TFT-LCD Touch Screen



iM8

10.4" Color TFT-LCD Touch Screen



No-fan Design



Pitch Tone



Thermal Recorder



Shortcut Menu

End Tidal Carbon Dioxide Monitoring



EDAN G2 CO₂ (sidestream)

Superior water trap design for accurate monitoring

iCARBTM algorithm with intelligent CO₂ pseudo wave identification technology

Multiple sampling accessories as options

Respironics CO₂ (mainstream/sidestream)

Plug & play module design

Dehumidification tube instead of water trap

Low sampling rate of 50ml/min

MASIMO Anesthetic Gas/O₂ (Mainstream/Sidestream)



Unique plug & play module for mainstream anesthetic gas monitoring

The no moisture sampling line design suitable for all clinical application

A stable low flow of 50ml/min


Multiple Display Modes



Standard Display



Large Font

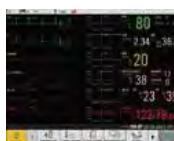

Trend View



OxyCRG



Back to Bed



Wall Mount



Rolling Stand